

See discussions, stats, and author profiles for this publication at: <https://www.researchgate.net/publication/262229214>

Rapid identification of bacteria in blood cultures by mass-spectrometric analysis of volatiles

ARTICLE *in* JOURNAL OF CLINICAL PATHOLOGY · MAY 2014

Impact Factor: 2.92 · DOI: 10.1136/jclinpath-2014-202301 · Source: PubMed

READS

22

7 AUTHORS, INCLUDING:



Simone Cristoni

ISB srl

52 PUBLICATIONS 556 CITATIONS

SEE PROFILE



Paolo Brambilla

Università degli Studi di Milano-Bicocca

115 PUBLICATIONS 3,200 CITATIONS

SEE PROFILE



Pablo Martinez-Lozano Sinues

ETH Zurich

39 PUBLICATIONS 452 CITATIONS

SEE PROFILE

Rapid identification of bacteria in blood cultures by mass-spectrometric analysis of volatiles

INTRODUCTION

Blood cultures are routine tests to determine whether micro-organisms have entered the patient's bloodstream. Automated systems, based on the detection of CO₂ increase in the culture media, have considerably improved the screening efficiency for the

Table 1 Volatile compounds detected in real-time in the head-space of the three bacteria investigated

m/z	Molecular formula	<i>Staphylococcus aureus</i> (n=15) mean (SD)	<i>Escherichia coli</i> (n=15) mean (SD)	<i>Streptococcus pneumoniae</i> (n=12) mean (SD)	Medium (n=15) mean (SD)	p-Value (ANOVA)
57.07028		1.46 (0.31)	1.45 (0.47)	1.89 (1.09)	1.11 (0.63)	1.9E-01
59.04946	C ₃ H ₇ O	6.03 (5.48)	6.57 (5.57)	3.41 (2.42)	3.81 (2.55)	2.2E-01
69.06999	C ₅ H ₉	5.53 (1.81)	7.47 (2.62)	8.82 (2.93)	5.48 (2.02)	4.9E-03
71.08584	C ₅ H ₁₁	0.94 (0.2)	2.32 (1.32)	1.19 (0.99)	0.51 (0.25)	7.7E-04
<i>73.04711</i>		<i>29.08 (15.25)</i>	<i>20.82 (15.1)</i>	<i>23.55 (12.98)</i>	<i>32.27 (10.25)</i>	<i>3.0E-01</i>
73.06511	C ₄ H ₉ O	16.05 (6.4)	30.47 (13.77)	30 (13.61)	13.79 (5.49)	2.2E-03
<i>89.04203</i>		<i>1.9 (1.64)</i>	<i>0.87 (0.99)</i>	<i>2.01 (1.45)</i>	<i>2.98 (2.61)</i>	<i>6.2E-02</i>
89.05991	C ₄ H ₉ O ₂	16.37 (7.53)	29.43 (23.98)	22.06 (6.88)	12.34 (18.6)	8.1E-02
<i>91.05756</i>	<i>C₄H₁₇S</i>	<i>57.46 (34.51)</i>	<i>40.69 (27.27)</i>	<i>32.04 (18.56)</i>	<i>67.07 (37.54)</i>	<i>6.7E-02</i>
95.0609		0.37 (0.33)	0.64 (0.32)	0.42 (0.26)	0.25 (0.33)	5.1E-02
95.08561		2.28 (0.96)	2.19 (1.03)	2.68 (0.76)	1.85 (0.74)	3.6E-01
97.06487	C ₆ H ₉ O	6.47 (1.39)	6.69 (1.64)	8.37 (3.42)	4.12 (1.46)	7.4E-02
105.0547	C ₄ H ₉ O ₃	4.13 (2.91)	9.17 (5.23)	10.85 (3.01)	0.43 (0.7)	1.5E-04
<i>105.073</i>	<i>C₅H₁₃S</i>	<i>61.17 (33.58)</i>	<i>28.99 (25.86)</i>	<i>63.24 (34.16)</i>	<i>73.26 (36.1)</i>	<i>8.2E-03</i>
105.0911	C ₅ H ₁₃ O ₂	0.89 (0.47)	1.42 (0.76)	1.22 (0.69)	0.5 (0.33)	9.3E-02
109.076	C ₆ H ₉ N ₂	2.9 (1.66)	3.86 (1.69)	3.74 (1.24)	0.95 (1.16)	2.1E-01
109.1012		4.41 (2.88)	10.47 (11.34)	11.17 (7.03)	2.53 (1.3)	5.5E-02
113.0961	C ₇ H ₁₃ O	8.78 (4.45)	16.48 (8.86)	18.26 (5.1)	3.99 (1.3)	9.2E-04
115.0754	C ₆ H ₁₁ O ₂	75.58 (23.7)	77.72 (18.78)	90.1 (13.65)	37.04 (20.75)	1.4E-01
119.0887	C ₆ H ₁₅ S	46.51 (21.06)	69.7 (32.42)	2.42 (2.92)	2.4 (1.52)	2.4E-08
119.1066		5.35 (2.25)	5.31 (1.72)	6.3 (1.99)	2.57 (1.26)	3.7E-01
127.1118	C ₈ H ₁₅ O	10.14 (3.71)	19.12 (14.96)	20.14 (7.1)	6.29 (2.21)	2.0E-02
131.0887	C ₇ H ₁₅ S	13.49 (8.18)	11.5 (7.63)	19.66 (20.32)	16.57 (10.68)	2.5E-01
142.0631		1.73 (2.29)	5.61 (5.31)	5.35 (4.24)	0.1 (0.39)	2.6E-02
<i>145.1043</i>	<i>C₈H₁₇S</i>	<i>3.84 (2.79)</i>	<i>3.6 (2.44)</i>	<i>3.14 (2.5)</i>	<i>5.51 (3.78)</i>	<i>7.8E-01</i>
145.1223	C ₈ H ₁₇ O ₂	9.92 (8.96)	28.41 (16.95)	28.85 (10.78)	3.42 (2.98)	2.7E-04
162.1489	C ₈ H ₂₀ NO ₂	3.33 (2.34)	5.87 (4.13)	5.46 (3.86)	2.92 (2.96)	1.2E-01
<i>163.0969</i>	<i>C₇H₁₅O₄</i>	<i>25.35 (8.43)</i>	<i>17.83 (8.83)</i>	<i>13.32 (5.03)</i>	<i>29.7 (12.1)</i>	<i>9.7E-04</i>
177.1485	C ₉ H ₂₁ O ₃	1 (1.35)	3.18 (2.29)	3.74 (2.9)	0.18 (0.29)	5.3E-03
203.1641	C ₁₁ H ₂₃ O ₃	0.99 (1.89)	2.3 (2.74)	2.42 (3.32)	0.71 (1.44)	2.9E-01
219.1745		0.51 (0.39)	1.41 (1.06)	1.76 (1.75)	0.31 (0.24)	1.8E-02
233.2112	C ₁₃ H ₂₉ O ₃	2.15 (3.2)	7.36 (6.24)	8.45 (7.07)	0.18 (0.36)	1.1E-02
308.2796	C ₁₄ H ₃₆ O ₃ N ₄	6.32 (12.01)	36.31 (39.34)	30.44 (41.73)	5.78 (12.84)	4.4E-02

Average signal of compounds in italics was found to be higher than for the three pathogens investigated.

and principal component analysis (PCA). A k-nearest-neighbour (k-nn, k=1) classifier was employed to predict sample class in a leave-one-out cross-validation (LOOCV).

RESULTS

An immediate mass spectral fingerprint was obtained from all samples. Thirty-three signals in the range 57–308 m/z were found to arise above the background upon injection of the gas. Table 1 lists their m/z values, molecular formula, mean (SD) and p value resulting from the ANOVA test. The average intensity of 27 out of the 33 compounds detected was higher in at least one of the bacterial strains than in the medium. Figure 1 shows an example of the real-time measurement of nine samples (three biological replicates of each bacteria strain) for six representative compounds. Note that each measurement requires a few seconds. Thus, in approximately 10 min, rich volatile fingerprints were obtained for the nine samples. Note how the signal rises upon the injection of the sample with a characteristic intensity, depending on the bacterial type, with a reasonable reproducibility for the three different biological replicates. For example, the compound at m/z 119.0887 (C₆H₁₅S) was present in *S. aureus* and *E. coli*, but not in *S. pneumoniae* (figure 1A) and the medium (table 1).

The corresponding PCA score plot (figure 2) of the mass spectral fingerprints reveals a clear clustering according to bacterial specimen.⁸ The LOOCV resulted in 40 out of 42 (95.2%) samples correctly classified. One *S. aureus* and one *S. pneumoniae* sample were misclassified as *E. coli*.

CONCLUSIONS

This exploratory study suggests that it is possible to identify bacteria by analysing the gases produced in the headspace of blood cultures by mass spectrometry in a rapid fashion with no sample preparation. Consistently with recent real-time mass spectrometric studies of bacterial cultures,^{3–4} relatively low molecular weight volatiles like acetone (m/z 59), acetoin (m/z 89) and butanal (m/z 73) were found in the cultures' headspace. In addition to this, SESI-MS captures rich fingerprints of volatiles expanding well beyond 300 Da.⁹ This data suggests that such high molecular weight vapours may be key to deliver high accuracy in pattern recognition of bacteria. For example, the ion detected at m/z 308.2796 (C₁₄H₃₆O₃N₄), was found to be the second largest relative contributor in PC3, which essentially separates *E. coli* from *S. aureus* and *S. pneumoniae* (see loadings plot in online supplementary figure S1). However, this pilot study has

detection of bacteria.¹ However, further identification of bacteria still requires time-consuming culturing procedures.

It has been suggested that along with CO₂, bacterial cultures emit characteristic volatile organic compounds that may be valuable for characterisation.² Recently, a number of technological developments have allowed the detection of trace gases with minimal or no sample preconcentration steps.^{3–5} In the current study, we investigated whether the volatiles emitted from bacterial cultures (mimicking routine clinical blood cultures) could be distinguished from each other in a rapid fashion by secondary electrospray ionisation-mass spectrometry (SESI-MS).^{6–7}

METHODS

10³ colony forming units were inoculated in 42 aerobic bottles: 15 with

Staphylococcus aureus ATCC 29213 strain, 15 with *Escherichia coli* ATCC 25922 and 12 *Streptococcus pneumoniae* ATCC 49619. They were incubated in the automatic instrument BacT/ALERT 3D (Biomérieux Clinical Diagnostic, France). The system automatically detects the positivity of the sample when the amount of CO₂ produced by bacterial growth reaches a given threshold. The bacterial identity was confirmed by culturing according to UK guidelines.

The headspace gas of the positive bottles was sampled with a 5 mL syringe. The gas was subsequently injected into a modified high-resolution mass spectrometer (Orbitrap, Thermo) located in the hospital premises. As a result, an immediate mass spectral fingerprint was obtained. The mass spectra were subjected to one-way analysis of variance (ANOVA) test

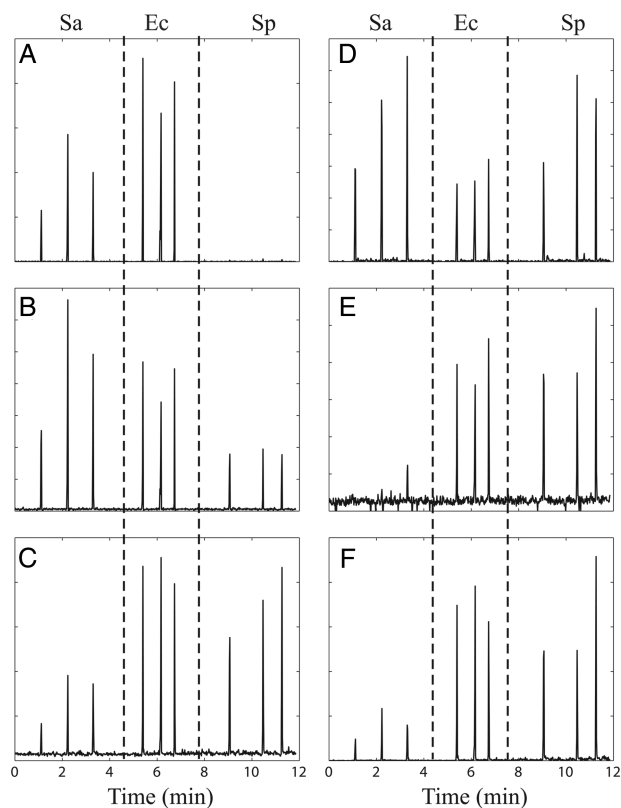


Figure 1 The headspace of commercial blood cultures can be sampled and mass analysed in a few seconds, providing an indication of the pathogen present. The plot displays the intensity as a function of time of six exemplary compounds detected in the headspace of the bacterial cultures. Three biological replicates of *Staphylococcus aureus* (Sa), *Escherichia coli* (Ec) and *Streptococcus pneumoniae* (Sp) were injected at minutes 1.1, 2.2, 3.3; 5.4, 6.1, 6.7 and 9.1, 10.5, 11.3, respectively (A) m/z 119.0887 ($C_6H_{15}S$); (B) m/z 163.0969 ($C_7H_{15}O_4$); (C) m/z 73.06511 (C_4H_9O); (D) m/z 105.073 ($C_5H_{13}S$); (E) m/z 219.1745; (F) m/z 308.2796 ($C_{14}H_{36}O_3N_4$).

some limitations. For example, during this exploratory phase, we wanted to determine the feasibility of this approach by using pure bacterial standards with relatively high colony-forming units. Further research is needed to confirm these findings in expectedly more complex matrices like blood cultures from patients.

Claudia Ballabio,¹ Simone Cristoni,² Giovanni Puccio,³ Malcolm Kohler,⁴ Maria Roberta Sala,¹ Paolo Brambilla,^{1,5} Pablo Martinez-Lozano Sinues^{6,*}

¹Clinical Pathology Department, Desio Hospital, Desio, Italy

²ISB Ion Source & Biotechnologies, Milan, Italy

³Emmanuele Scientific Research Association, Palermo, Italy

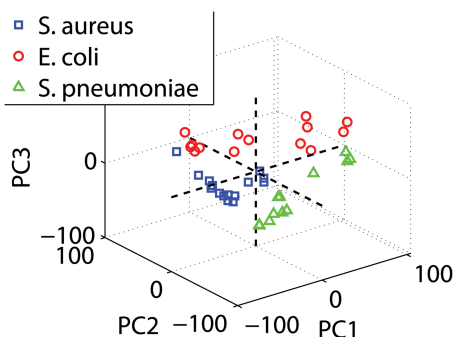


Figure 2 Proof-of-principle of rapid characterisation of three bacterial strains grown in blood bacterial culture vials. Score plot of principal component analysis of the mass spectral fingerprint of the volatiles emitted by three different bacterial agents: *Staphylococcus aureus* (n=15 samples), *Escherichia coli* (n=15) and *Streptococcus pneumoniae* (n=12). Clustering according to pathogen is evident, suggesting that before disposal of the positive vials, the analysis of its headspace may be a valuable source of information. Note that each of the measurements shown in the graph requires just a few seconds.

⁴Pulmonary Division, University Hospital Zurich, Zurich, Switzerland

⁵Experimental Medicine Department, School of Medicine, University of Milano-Bicocca, Monza, Italy

⁶National Research Council-Institute for Biomedical Technologies, Segrate, Italy

*Current address: Department of Chemistry and Applied Biosciences, ETH Zurich, Zurich, Switzerland

Correspondence to Dr Pablo Martinez-Lozano Sinues, ETH Zurich, Department of Chemistry and Applied Biosciences, Vladimir-Prelog-Weg 3, Zurich CH-8093, Switzerland; pablo.mlsinues@org.chem.ethz.ch

Contributors SC, PB and PM-LS designed the study. CB, SC, MRS and PM-LS conducted the experiments. PM-LS analysed the data. SC, GP, MK, MRS, PB and PM-LS wrote the manuscript.

Funding This research was supported by a Marie Curie Grant (PMLS) within the 7th European Community Framework Programme (PIEF-GA-2008-220511).

Competing interests None.

Provenance and peer review Not commissioned; externally peer reviewed.

► Additional material is published online only. To view please visit the journal online (<http://dx.doi.org/10.1136/jclinpath-2014-202301>).



CrossMark

To cite Ballabio C, Cristoni S, Puccio G, et al. *J Clin Pathol* 2014;**67**:743–746.

Received 17 March 2014

Revised 11 April 2014

Accepted 16 April 2014

Published Online First 9 May 2014

J Clin Pathol 2014;**67**:743–746.

doi:10.1136/jclinpath-2014-202301

REFERENCES

- 1 Thorpe TC, Wilson ML, Turner JE, et al. BacT/Alert: an automated colorimetric microbial detection system. *J Clin Microbiol* 1990;**28**:1608–12.
- 2 Bos LDI, Sterk PJ, Schultz MJ. Volatile metabolites of pathogens: a systematic review. *PLoS Pathog* 2013;**9**: e1003311.
- 3 Bunge M, Araghipour N, Mikoviny T, et al. On-line monitoring of microbial volatile metabolites by proton transfer reaction-mass spectrometry. *Appl Environ Microbiol* 2008;**74**:2179–86.
- 4 Sovova K, Cepl J, Markos A, et al. Real time monitoring of population dynamics in concurrent bacterial growth using SIFT-MS quantification of volatile metabolites. *Analyst* 2013;**138**: 4795–801.
- 5 Chippendale TWE, Gilchrist FJ, Spanel P, et al. Quantification by SIFT-MS of volatile compounds emitted by in vitro cultures of *S aureus*, *S pneumoniae* and *H influenzae* isolated from patients with respiratory diseases. *Anal Methods* 2014;**6**:2460–72.
- 6 Dillon LA, Stone VN, Croasdell LA, et al. Optimisation of secondary electrospray ionisation (SESI) for the trace determination of gas-phase volatile organic compounds. *Analyst* 2010;**135**:306–14.
- 7 Martinez-Lozano P, Rus J, Fernández de la Mora G, et al. Secondary Electrospray Ionization (SESI) of ambient vapors for explosive detection at concentrations below parts per trillion. *J Am Soc Mass Spectrom* 2009;**20**:287–94.

PostScript

- 8 Zhu J, Bean HD, Kuo Y-M, *et al.* Fast detection of volatile organic compounds from bacterial cultures by secondary electrospray ionization-mass spectrometry. *J Clin Microbiol* 2010;48:4426–31.
- 9 Martínez-Lozano P, Fernández de la Mora J. Electrospray ionization of volatiles in breath. *Int J Mass spectrom* 2007;265:68–72.



Rapid identification of bacteria in blood cultures by mass-spectrometric analysis of volatiles

Claudia Ballabio, Simone Cristoni, Giovanni Puccio, et al.

J Clin Pathol 2014 67: 743-746 originally published online May 9, 2014
doi: 10.1136/jclinpath-2014-202301

Updated information and services can be found at:
<http://jcp.bmj.com/content/67/8/743.full.html>

These include:

Data Supplement

"Supplementary Data"

<http://jcp.bmj.com/content/suppl/2014/05/09/jclinpath-2014-202301.DC1.html>

References

This article cites 9 articles, 3 of which can be accessed free at:

<http://jcp.bmj.com/content/67/8/743.full.html#ref-list-1>

Email alerting service

Receive free email alerts when new articles cite this article. Sign up in the box at the top right corner of the online article.

Notes

To request permissions go to:

<http://group.bmj.com/group/rights-licensing/permissions>

To order reprints go to:

<http://journals.bmj.com/cgi/reprintform>

To subscribe to BMJ go to:

<http://group.bmj.com/subscribe/>

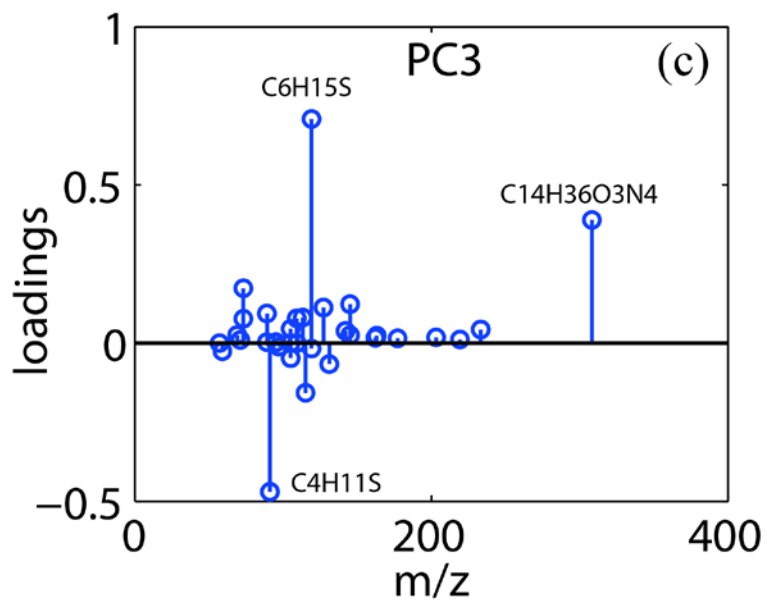
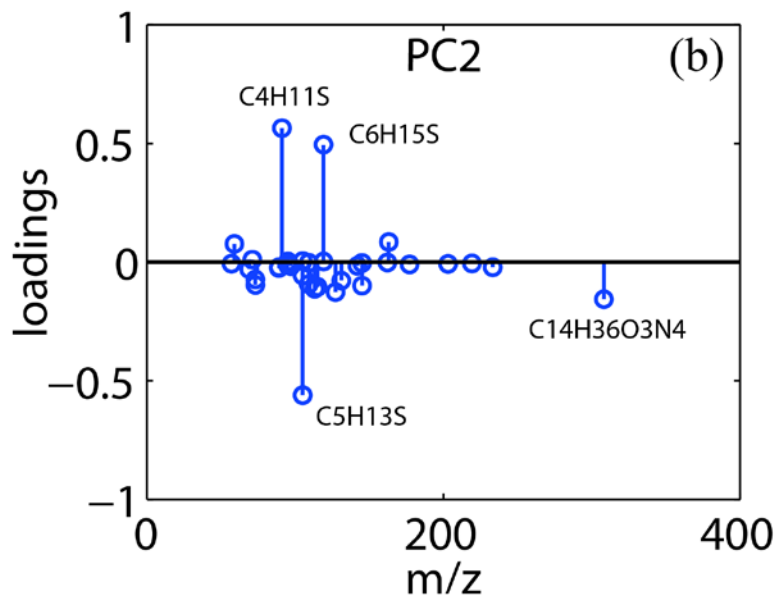
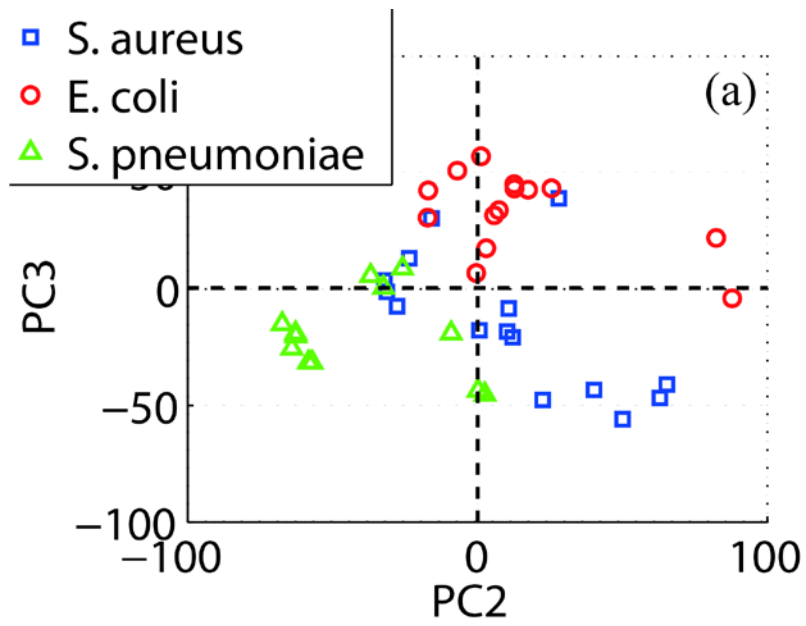


Figure S1. (a) Score plot of PC2 vs. PC3 suggests a clustering according to bacterial strain. The corresponding loadings plots of PC2 (b) and PC3 (c).